




HIGH PERFORMANCE DIGITAL IMAGING
made easy

MICROPUBLISHER 5.0 & 3.3 RTV – COOLED OR NON-COOLED

High-Resolution IEEE 1394 FireWire™ Digital CCD Color Camera with High-Speed Real-Time Viewing

The QImaging MicroPublisher with Real-Time Viewing (RTV) delivers unsurpassed interactivity and productivity by combining ultra-high-resolution images with video-like, full-field-of-view frame rates up to 30fps. Scanning, framing, and focusing have never been easier than with the MicroPublisher RTV. The 30-bit color digitization produces high-quality images of brightfield, darkfield, and fluorescence work. For demanding low-light applications, the MicroPublisher RTV Cooled camera minimizes thermal noise during long exposure times. With an IEEE 1394 FireWire™ digital interface, the MicroPublisher RTV is easy to install, requiring a single wire to connect the camera to a computer or laptop. The MicroPublisher RTV eliminates expenses, installation problems, and inconveniences associated with framegrabbers and external power supplies. All cameras ship with image-capture software. A large selection of specialty software applications is available from QImaging's software partners. A **Software Development Kit (SDK)** is available upon request for interfacing the MicroPublisher RTV with custom software.

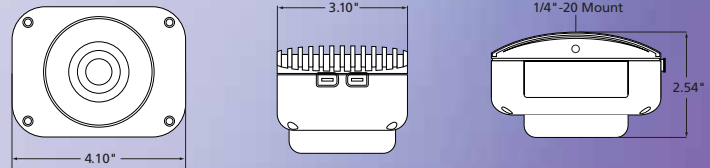

Note: Microscope is shown for illustration only and is not included.

CAMERA MODELS

Includes: IEEE 1394 FireWire™ cable, IEEE 1394 PCI card, QCapture software, QCapture Pro software, & access to SDK

MicroPublisher 5.0 RTV Cooled

Model: MP5.0-RTV-CLR-10-C
CCD Digital Camera, Color,
30 Bits with Peltier Cooling

MicroPublisher 5.0 RTV Non-Cooled

Model: MP5.0-RTV-CLR-10
CCD Digital Camera, Color, 30 Bits

MicroPublisher 3.3 RTV Cooled

Model: MP3.3-RTV-CLR-10-C

MicroPublisher 3.3 RTV Non-Cooled

Model: MP3.3-RTV-CLR-10

FEATURES

High-Resolution,
5- or 3.3-Million-Pixel Sensor

Real-Time Viewing (RTV)

Flexible Exposure Control
from 1.6ms to 17.9min

Peltier Cooling

ROI (Region of Interest)

Binning

IEEE 1394 FireWire™

Extensive Third-Party
Software Support

BENEFITS

- Highly detailed, sharp images suitable for publication

- Previewing & focusing in real time
- 30fps (full field of view) with MicroPublisher 3.3 RTV
- 25fps (full field of view) with MicroPublisher 5.0 RTV

- Optimal integration over a wide range of light levels

- Minimizes thermal noise during low-light imaging

- Higher frame rates for previewing & focusing

- Increases sensitivity for quantitation & imaging of very low light levels
- Increases frame rate

- Simple connectivity
- Ease of use & installation
- Portability with laptop computer
- Simultaneous use of multiple cameras through a single port
- Single-cable operation (no external power supply or control unit)

- Choose from a large selection of life science & industrial software for microscopy, machine vision, & video-streaming applications



MICROPUBLISHER 5.0 & 3.3 RTV SPECIFICATIONS

APPLICATIONS

High-resolution still images for publication, documentation, and archiving in:

- Brightfield, Phase-Contrast, & Darkfield Microscopy
- Fluorescence Microscopy
- Pathology
- Histology
- Cytology
- Hematology
- Document Imaging
- Still-Image Animation

CCD SENSOR

Light-Sensitive Pixels MicroPublisher 5.0 RTV MicroPublisher 3.3 RTV	5 million real pixels; 2560 x 1920 3.3 million real pixels; 2048 x 1536
Binning Modes	2x2, 3x3, 4x4 in full color
ROI (Region of Interest)	From 1x1 pixels up to full resolution, continuously variable in single-pixel increments
Exposure/Integration Control	1.6ms to 17.9min in 1µs increments
Sensor Type MicroPublisher 5.0 RTV MicroPublisher 3.3 RTV	Sony® ICX282 progressive-scan interline CCD (color) Sony® ICX252 progressive-scan interline CCD (color)
Pixel Size MicroPublisher 5.0 RTV MicroPublisher 3.3 RTV	3.4µm x 3.4µm 3.45µm x 3.45µm
Cooling Available	Yes (optional)
Cooling Type	Peltier thermoelectric cooling to 10°C below ambient
Digital Output	10 bits
Readout Frequency	20, 10, 5, 2.5MHz
Frame Rate MicroPublisher 3.3 RTV MicroPublisher 5.0 RTV	30fps full field of view (higher fps with ROI functions) 25fps full field of view (higher fps with ROI functions)

CAMERA

Computer Platforms/Operating Systems	Windows® & Mac OS*
Digital Interface	IEEE 1394 FireWire™
Shutter Control	Electronic shutter, no moving parts
Trigger Types	Internal, Software
Optical Interface MicroPublisher 5.0 RTV MicroPublisher 3.3 RTV	2/3", C-mount optical format 1/2", C-mount optical format
Threadmount	1/4" — 20 mount
Power Requirements	3.8W (non-cooled); 6.7W (cooled); 8-24V
Weight	710g
Warranty	2 years
Operating Environment	0 to 35°C (32 to 95°F)
Humidity	Less than 80% non-condensing at 35°C (95°F)

*Refer to QImaging website for detailed listing of supported operating systems.
Note: Specifications are nominal and subject to change.

ISO 9001:2000



04-0015A-C



Tel 604.708.5061
Fax 604.708.5081
INFO@QIMAGING.COM
WWW.QIMAGING.COM

FireWire and Mac OS are trademarks of Apple Computer, Inc., registered in the U.S. and other countries. Sony is a registered trademark of Sony Corporation. Windows is a registered trademark of Microsoft Corporation in the United States and other countries. Other brand and product names are the trademarks or registered trademarks of their respective owners and manufacturers.